


RETI-CRYO

GENERAL INFORMATION

Reti-Cryo is an electrically controlled, easy-to-use cryo surgical unit for ophthalmic surgical needs. The foot pedal is pressed for the desired freeze period and the system automatically initiates the de-freeze cycle. To meet various defreeze cycle durations, a front panel setting allows to set the defreeze duration. These units are supplied with retinal probes or optional probe for glaucoma etc.

FEATURES & SPECIFICATIONS

- Quick Freeze and De-freeze cycles
- Easy to operate electrical foot switch
- Electrically controlled N₂O Cryo Surgical unit
- Automatic De-Freeze cycle initiation when the foot pedal is released
- Adjustable De-Freeze timing with control switches in the front panel
- Electrical power: AC 100-240 Volts, 50/60 Hz, Single Phase
- Weight: 4 Kg / Dimensions: L 290 mm x D 256 mm x H 62 mm

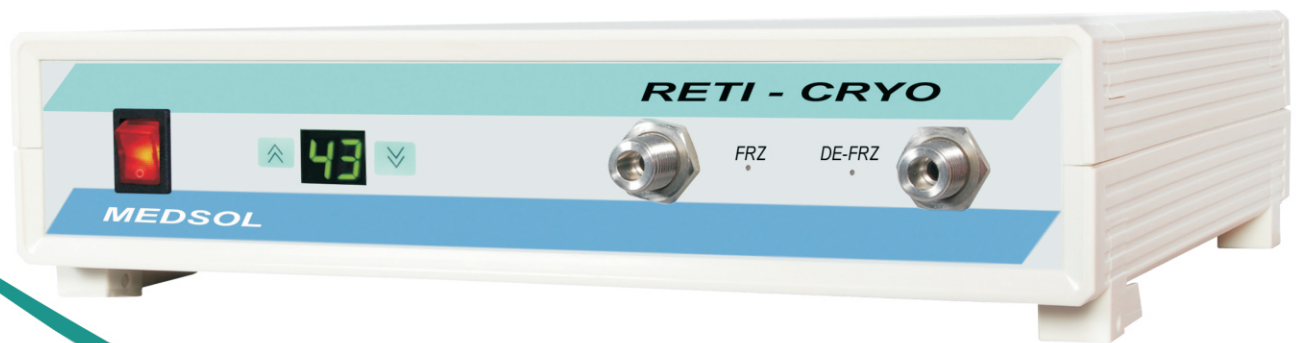

Manufactured by:



ICON MEDICARE

#304-A/2, TNHB, Phase VII, Brindavan Nagar,
HOSUR - 635 109. Tamilnadu, INDIA.

Phone : +91 4344 243313

E-mail : info@iconmedicare.com / www.iconmedicare.com

NEW DELHI
Mob.: +91 9944678009

MUMBAI
Mob.: +91 9323284971

BANGALORE
Mob.: +91 8220009934

CHENNAI
Mob.: +91 8144433313

KOLKATTA
Mob.: +91 9842191245

AHMEDABAD - Mob.: +91 9327065024

HYDERABAD - Mob.: +91 9000929890

BI-POLAR COAGULATOR /CAUTERY

SURGICAL EQUIPMENT

WET-FIELD BI-POLAR

GENERAL INFORMATION

Bi-Polar Cautery is a compact and portable unit that is designed for use in cauterizing the tissues during the ophthalmic surgery. This Bipolar cautery allows fixed power delivery and the user can set the necessary power using the front panel control.

FEATURES & SPECIFICATIONS

- Fixed power delivery
- Output short circuit protection
- Compact, lightweight and portable
- 10 Watts high precision power output
- Autoclavable Bi-polar Forceps or Pencil Eraser
- Variable power setting from 05-100% at a step of 05%
- Electrical power: AC100-240V, 50/60 Hz, Single Phase
- Weight: 1Kg / Dimensions: L 160 mm x D 197 mm x H 50 mm



Manufactured by:



ICON MEDICARE

#304-A/2, TNHB, Phase VII, Brindavan Nagar,
HOSUR - 635 109. Tamilnadu, INDIA.

Phone : +91 4344 243313

E-mail : info@iconmedicare.com / www.iconmedicare.com

NEW DELHI
Mob.: +91 9944678009

MUMBAI
Mob.: +91 9323284971

BANGALORE
Mob.: +91 8220009934

CHENNAI
Mob.: +91 8144433313

KOLKATTA
Mob.: +91 9842191245

AHMEDABAD - Mob.: +91 9327065024

HYDERABAD - Mob.: +91 9000929890